

SINGLE POINT STIRRERS WITH INTERNAL CONTROL


Designed for chemical, biotechnical and medical applications of all kinds, developed for the daily lab use and continuous operation.

MIX 1

Ultra flat, wear-free magnetic stirrer with 1 stirring position, up to 10,000 ml stirring volume, also suitable for viscous media.

- ⌘ Maintenance- and wear-free
- ⌘ Up to 10 litres or for viscous media
- ⌘ Outstanding long-life quality
- ⌘ Developed and made in Germany

magnetic motion

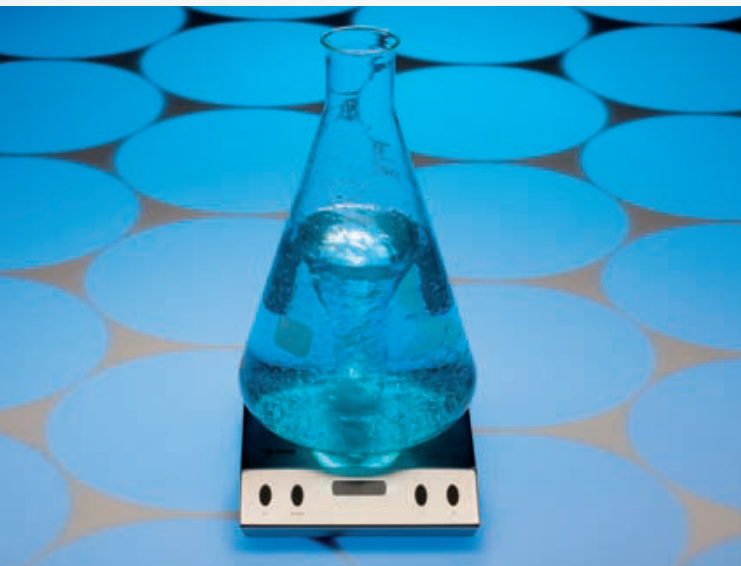

Type	MIX 1
Stirring points	1
Stirring volume/stirring point	1 - 10,000 ml
Speed range	100 - 2,000 rpm
Power setting	5/10/15/20 W (4-step)
Material housing	Stainless steel
Permitted operation conditions	-10 °C up to +50 °C (at 80% humidity)
Electrical data	100-240 V / 50-60 Hz / 1.5 A
Dimensions (WxDxH)	188 x 238 x 38 mm
Protection category	IP64
Weight (gross)	approx. 4 kg
Order no.	30001

3 years warranty on material and manufacturing.
We reserve the right to make technical alterations. We do not assume liability for errors in printing or typing.

100% maintenance-free and wear-free inductive 2mag-Magnetic-Drive concept for powerful mixing of solutions, speed range from 100 up to 2,000 rpm, jerk-free stirring even at low speeds.

4-step power setting for high power with regard to large and viscous stirring volumes and reduced power for the warming-free continuous operation.

Clear digital display for speed- and power settings, SoftStart for reliable catching/centering and safe accelerating of the stirring bar. QuickSet for quick setting of start- and maximum speed, automatic memory function to save the last settings, outstanding long-life cycle of the hermetically sealed drive system, resistant stainless steel housing for easy and fast cleaning.

Accessories

ASTEROID 40 - Stirring bar
Stirring bar remover PTFE

