

VERTICAL LAMINAR FLOW FOR SAMPLE PROTECTION

# CRUMA670FL



The CRUMA vertical laminar flow hood with ISO Class 5 air purity creates a sterile environment free of particles and/or microorganisms that ensures the reliability of any test inside the hood. It incorporates a HEPA H-14 absolute filter whose minimum efficiency is of 99.995% for particles of 0.3  $\mu\text{m}$ , and a 15W UV germicide lamp designed to decontaminate the hood before and after each use.

It includes **new important technical features** in order to satisfy the operator in his everyday task and his protection at work, taking care of the environment as well. All this with the same quality as always and up to 7 years guarantee.



## TECHNICAL FEATURES

Number of filters	1
Number of IP44 fans	1
Average face velocity	0,40 m/s
Total electrical power consumption	65 W
Voltage-Frequency	110-220 V / 50-60 Hz
LED light intensity	1500 lm
UV germicidal lamp	15W
Noise level	59 dB
<b>Packaging:</b> 100% recycled wooden box with international phytosanitary certificate	Volume 0,28 m <sup>3</sup> Weight 65 Kg

## SIZES (MM)


External			Internal		
Width	Depth	Height	Width	Depth	Height
600	600	937	576	556	620

Filters with microchip and USB connection, *intelligent filtering*

Information storage and interaction with the alarm system.



## SIZES (mm)

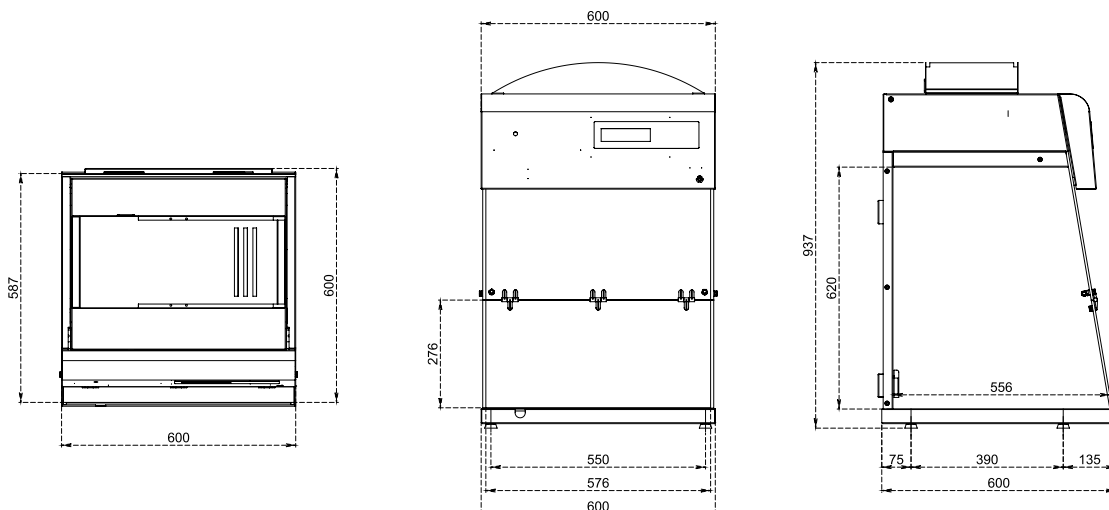

VERTICAL LAMINAR FLOW FOR SAMPLE PROTECTION

# CRUMA870FL



The **CRUMA vertical laminar flow hood** with ISO Class 5 air purity creates a sterile environment free of particles and/or microorganisms that ensures the reliability of any test inside the hood. It incorporates a HEPA H-14 absolute filter whose minimum efficiency is of 99.995% for particles of 0.3  $\mu\text{m}$ , and a 15W UV germicide lamp designed to decontaminate the hood before and after each use.

It includes **new important technical features** in order to satisfy the operator in his everyday task and his protection at work, taking care of the environment as well. All this with the same quality as always and up to 7 years guarantee.




## TECHNICAL FEATURES

Number of filters	1
Number of IP44 fans	1
Average face velocity	0,45 m/s
Total electrical power consumption	100 W
Voltage-Frequency	110-220 V / 50-60 Hz
LED light intensity	1500 lm
UV germicidal lamp	15W
Noise level	57 dB
<b>Packaging:</b> 100% recycled wooden box with international phytosanitary certificate	Volume 0,56 m <sup>3</sup> Weight 92 Kg

## SIZES (MM)

External			Internal		
Width	Depth	Height	Width	Depth	Height
800	600	1178	775	556	720

Do you have any *questions* or need *technical assistance*?

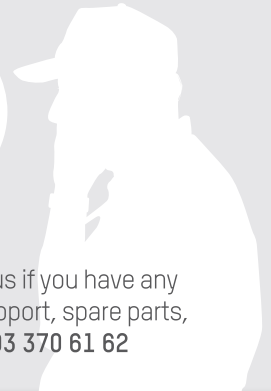

Contact your distributor or call us if you have any questions or need technical support, spare parts, maintenance service... ☎ +34 93 370 61 62

## SIZES (mm)

